

DUAL CLICK-TITE™ CONTAINERS

Specimen Collection and Transport System

Therapak's Dual Click-Tite™ Containers are ideal for maintaining the integrity of biological specimens during transit and storage.

Features and Benefits:

Durable Design

Composed of 100% virgin FDA approved high clarity polypropylene, these sturdy yet lightweight containers are 95kPa certified for air transport.

Leak-Resistant Dual Click-Tite™ Cap Seal

All containers feature fast screw threads with low-friction sealing surfaces for quick and easy container open and closure. Dual Click-Tite caps “click” twice when completely closed, indicating the container has been closed securely and is now leak-resistant. Inadvertent back-off of the cap during handling and transport is prevented by an anti-back-off system designed into the containers.

Freight Efficient Packaging


The 60mL x 53mm Dual Click-Tite containers are packed 560 per case in the same carton that holds 400 of the 90mL x 53mm containers. Choosing a 60mL x 53mm for your application not only reduces freight cost per container in this comparison – reduced container weight resulting from collecting and transporting less sample may also prove beneficial for inbound freight savings.

Fast, Accurate Temperature Reading

Containers are of a robust design while utilizing the absolute minimum of materials and energy during their production. Thinner sidewalls permit rapid temperature transfer for applications requiring immediate and accurate temperature sensing of the specimen.

Sterile and Tamper Evident

Sterile containers are complete with tamper evident label. Containers are stackable, graduated and include top and bottom knurls for ease of gripping with gloves.

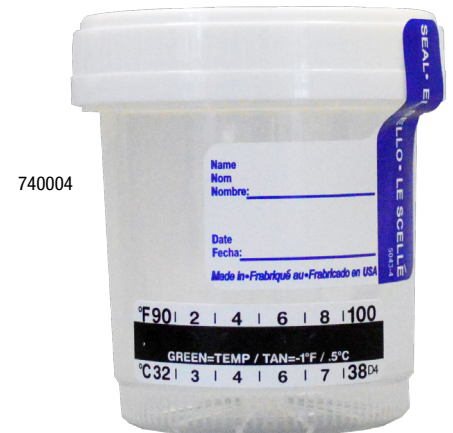

| Part No. | Item | Packaging | Inner Pkg |
|----------|--|-----------|-----------|
| 740001 | 60mL Sterile White Cap w/ Label | 560/cs | 4 x 140 |
| 740002 | 60mL Sterile White Cap w/ Label and Temp | 560/cs | 4 x 140 |
| 740003 | 90mL Sterile White Cap w/ Label | 400/cs | 4 x 100 |
| 740004 | 90mL Sterile White Cap w/ Label and Temp | 400/cs | 4 x 100 |
| 740005 | 90mL Sterile White Cap w/ Label and Temp | 400/cs | 8 x 50 |
| 740006 | 120mL Sterile White Cap w/ Label | 300/cs | 4 x 75 |

All containers are 53mm in diameter.

FOR MORE INFORMATION

(909) 267-2000/(770) 614-2931

info@therapak.com

www.therapak.com