

## Dual Click-Tite™ Containers

Therapak’s Dual Click-Tite™ Containers are ideal for maintaining the integrity of biological specimens during transit and storage.

### Durable Design

Composed of 100% virgin FDA approved high clarity polypropylene, these sturdy yet lightweight containers are 95kPa certified for air transport.

### Leak-Resistant Dual Click-Tite™ Cap Seal


All containers feature fast screw threads with low-friction sealing surfaces for quick and easy container open and closure. Dual Click-Tite caps “click” twice when completely closed, indicating the container has been closed securely and is now leak-resistant. Inadvertent back-off of the cap during handling and transport is prevented by an anti-back-off system designed into the containers.

### Fast, Accurate Temperature Reading

Containers are of a robust design while utilizing the absolute minimum of materials and energy during their production. Thinner sidewalls permit rapid temperature transfer for applications requiring immediate and accurate temperature sensing of the specimen.

### Sterile and Tamper Evident

Sterile containers are complete with tamper evident label. Containers are stackable, graduated and include top and bottom knurls for ease of gripping with gloves.

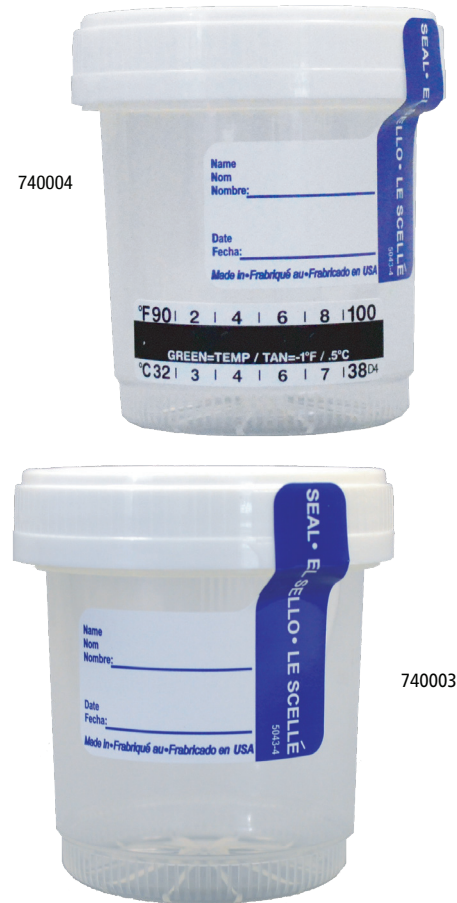

Part No.	Item	Packaging	Inner Pkg
740085	60mL Sterile White Cap with Label	500/cs	5 x 100
740086	60mL Sterile White Cap with Label and Temperature Strip	500/cs	5 x 100
740003	90mL Sterile White Cap with Label	400/cs	4 x 100
740004	90mL Sterile White Cap with Label and Temperature Strip	400/cs	4 x 100
740005	90mL Sterile White Cap with Label and Temperature Strip	400/cs	8 x 50
740006	120mL Sterile White Cap with Label	300/cs	4 x 75
740095	120mL Sterile White Cap with Label	300/cs	4 x 75

### FOR MORE INFORMATION

(909) 267-2000/(770) 614-2931

info@therapak.com

www.therapak.com