

USB Temperature Data Logger

Microstar™ USB Temperature Data Loggers provide an inclusive solution for monitoring temperature sensitive products. This single use device starts recording 10 minutes after activation, and will record more than 12,000 temperature data points* over the next 90 days unless manually stopped.


Activate the logger inside the waterproof packaging. Use the attached adhesive strip to affix the logger inside the container, or simply place it with your payload. To access the data, remove the packaging and insert into the USB slot of your computer. A fully detailed PDF report is generated automatically and includes specific logger information, temperature graph and data points.

Each logger is tested and validated by an accredited laboratory and traceable to NIST.

Features and Benefits:

- Single-use and easy to activate
- No additional software/hardware required[†]
- Temperature Range: -30°C to +70°C (-22°F to +158°F)
- Full data report provided including graph and data points

*All logger data is in Universal Coordinated Time (UTC)
[†]Adobe Reader is required to view report. Free download at adobe.com

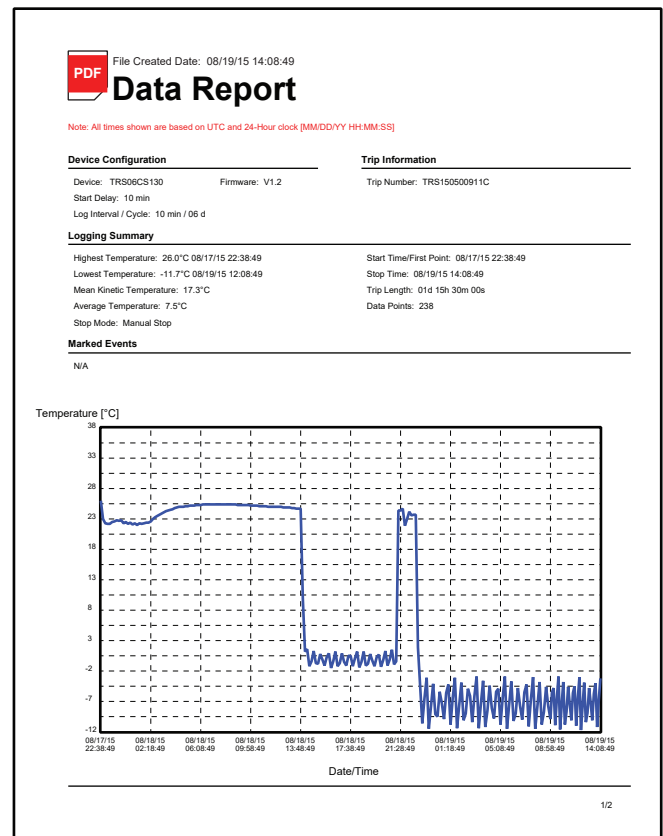

Part No.	Item	Temperature Range	Reading Interval	Maximum Duration	Packaging
365161	Microstar™ USB Temperature Data Logger	-30°C to +70°C	10 minutes	90 days	Each

Microstar™ is a trademark of Therapak, LLC.

FOR MORE INFORMATION

(909) 267-2000/(770) 614-2931

info@therapak.com

www.therapak.com