


DURAMARK[™] VIP SHIPPERS

Featuring NanoCool[®] technology

- **Effective:** Temperature tested to ISTA 7D standards
- **Efficient:** Uses less dry ice than conventional shippers
- **Economical:** Reduces freight costs with smaller shipper size (low DIM weight)
- **Environmentally friendly:** Packaging requires less storage space and creates less waste

The Duramark[™] Vacuum Insulated Panel (VIP) Shipper is ideal for the transport of frozen (<-20°C) shipments.

This sturdy shipper has been drop tested under dry ice conditions and temperature tested to ISTA 7D standards. The innovative insulation panels require less dry ice while offering long lasting temperature control. There is a range of options to meet your needs: the 1.5L VIP shipper can maintain temperatures below -20°C for 2 days using only 3 pounds of dry ice while the 10.4L size maintains -20°C for 6 days with 13 pounds of dry ice. The compact packaging requires less storage space, decreases freight costs and generates less waste.

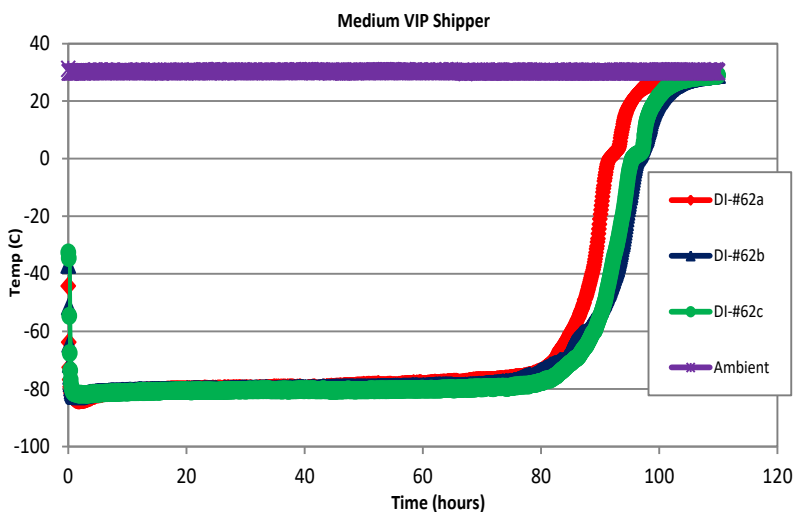

3.2L VIP Shipper product performance with 4.7 pounds of Dry Ice at a constant 30°C ambient temperature.

Part No.	Item	Duration below -20°C (hours) / Dry ice used (pounds)	Inner Dimensions	Outer Dimensions
562435	VIP Shipper - 1.5 Liter Capacity	48 / 2.8	8.4 x 4.5 x 2.5"	11.8 x 8.2 x 5.4"
562436	VIP Shipper - 3.2 Liter Capacity	72 / 4.3	8.4 x 4.4 x 5.3"	11.3 x 7.8 x 7.5"
562437	VIP Shipper - 8.2 Liter Capacity	144 / 11.0	9.1 x 8.6 x 6.4"	12.3 x 11.5 x 10"
562438	VIP Shipper - 10.4 Liter Capacity	144 / 13.0	9.1 x 8.6 x 8.1"	12.3 x 11.5 x 11.7"
562907	VIP Shipper - 14.7 Liter Capacity		11.3 x 8.9 x 8.1"	14.0 x 12.9 x 12.2"
562908	VIP Shipper - 18.5 Liter Capacity		11.3 x 8.9 x 10.4"	14.0 x 12.9" x 14.3"

FOR MORE INFORMATION

(909) 267-2000 / (770) 614-2931

info@therapak.com

www.therapak.com